

ACTIVATE YOUR ENDODONTICS

This innovative device agitates irrigation solutions during endodontic treatment. Evidence-based endodontics has shown that cavitation and acoustic streaming significantly improve debridement and the disruption of the smear layer and biofilm. The EndoActivator[®] System is designed to safely and vigorously energize the hydrodynamic phenomenon. Activated fluids promote deep cleaning and disinfection into lateral canals, fins, webs, and anastomoses. A cleaned root canal system facilitates 3-D obturation and long-term success. Activate your endodontic solutions today with this easy-to-use technology!

DESIGNED BY
Clifford J. Ruddle DDS
Robert H. Sharp DDS
Pierre Machtou DDS



EndoActivator Tips

CLINICAL FEATURES

- Safe and efficient technology
- Ergonomic handpiece design
- Strong, flexible polymer tips
- Simple, affordable system



EndoActivator
Handpiece/Driver


ENDO ACTIVATOR[®]

SAFE • EFFICIENT • SIMPLE

SYSTEM SIMPLICITY

The EndoActivator Tips are 22 mm in length and are available in three sizes: small, medium, and large. Specific features include:

- Snap-on/snap-off design
- Strong, flexible, medical-grade polymer composition
- Color-coded by size for identification
- Convenient depth gauge rings at 18, 19 and 20 mm
- Sold separately in reorder packs or available as a complete system with the Handpiece/Driver

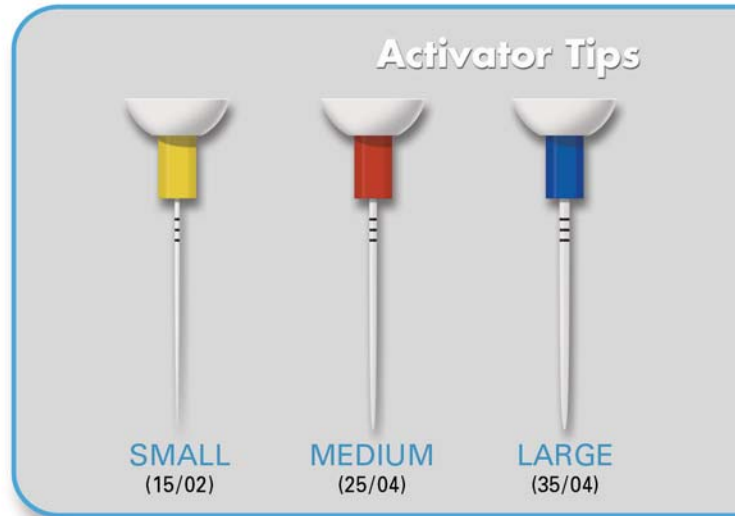

Handpiece/Driver



Rx Only

Artwork and images courtesy of Dr. Cliff Ruddle and Advanced Endodontics. EndoActivator[®] is a registered trademark of Endo Inventions LLC. SSEA Rev. 3 6/12

1-800-662-1202
For the latest information consult
www.TulsaDentalSpecialties.com

DENSPLY
TULSA DENTAL
SPECIALTIES